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• MAINTENANCE

Lubricate periodically with a drop of oil (Teflon) inside the connector for the plastic pipe in order to improve the movement of the knob, if necessary.

Characteristics:

- temperature of use: -20°C÷60
- storage temperature: -40°C ÷ 70°C
- level of protection of the manometer: IP 55
- scale division of the manometer: 50hPA/
- Precision class: 2,5

GUARANTEE

PROSIMED guarantees that the product LCLC1330CE is free from manufacturing or materials defects for a period of 12 months after the purchase. The manufacturer declines any guarantee for operations, such as cleaning, use and conservation. This guarantee is not valid under the following conditions:

- if the product is damaged for improper use or negligence;
- if the product is modified or repaired by unauthorised people;
- if the identification numbers of the product are illegible or altered.

RESPONSIBILITY

The manufacturer or retailer decline any responsibility for any loss and direct or indirect damages, expenses that are directly or indirectly connected to the use of the appliance LCLC1330CE or because of consequences of events, such as diagnosis, medical procedures and treatments taken by the medical personnel that may cause damages to persons or objects, or cause expenses. Inside the packaging, we supply the form “Warning of Malfunctioning” that must be filled and sent to the manufacturer for any accident or near miss that may occur during the intervention.

The code LCLC1330CE is supplied with the connections ① and ② assembled and it is ready to be applied upon carriages. In case of use of the wall regulator, dismantle the connection ① (key 17) and assembly it again in place of the connection ② (key 14) and, in place of the connection ① assembly the desired wall connection ③.

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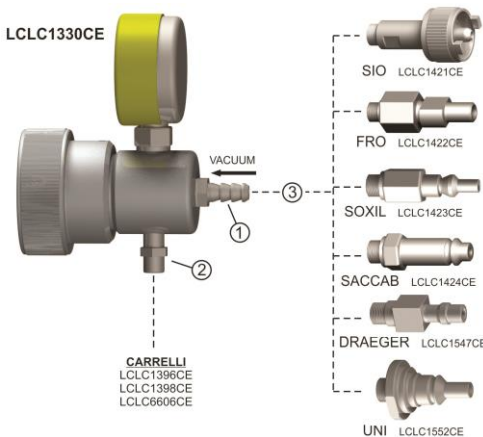
(EN) INSTRUCTIONS OF USE

THIS DEVICE MUST BE USED BY QUALIFIED PERSONNEL WHICH KNOW THE TECHNICAL CHARACTERISTICS OF THE VACUUM PLANT.

• **DESCRIPTION** The appliance LCLC1330CE is a regulator/gauge with vacuum measurement inside the plant.

• **INDICATIONS** The appliance LCLC1330CE connects the vacuum from a wall outlet to the containers of pockets for sucked liquids, having another possibility of regulation and measurement of vacuum inside the plant.

• OPERATING TEST TO BE IMPLEMENTED BEFORE EACH USE

1. Check that the wall outlet of connection is the one of vacuum (Yellow colour)
 2. Make sure that the source plant of vacuum is operating.
 3. Check that the Regulator is in position OFF (Close)
 4. Insert the regulator in the wall outlet
 5. Slowly turn to position ON
 6. Check that the needle of the device starts measuring
 7. Connect the connection pipe to the collection containers of sucked liquids
 8. Regulate the pressure of vacuum in a suitable position.
 9. Close the valve inside the collection containers of sucked liquids.
- Now the plant is ready to be used by medical personnel, which simply have to open the valve of containers to have the vacuum available.

• CLEANING

Use non-alcoholic detergent and dry well.
DO NOT USE CHLORINE-BASED PRODUCTS.
DO NOT IMMERSE

• MAINTENANCE

Lubricate periodically with a drop of oil (Teflon) inside the connector for the plastic pipe in order to improve the movement of the knob, if necessary.
Characteristics:
- temperature of use: -20°C÷60
- storage temperature: -40°C ÷ 70°C
- level of protection of the manometer: IP 55
- scale division of the manometer: 50hPa/
- Precision class: 2,5

GUARANTEE

PROSIMED guarantees that the product LCLC1330CE is free from manufacturing or materials defects for a period of 12 months after the purchase. The manufacturer declines any guarantee for operations, such as cleaning, use and conservation. This guarantee is not valid under the following conditions:
- if the product is damaged for improper use or negligence;
- if the product is modified or repaired by unauthorised people;
- if the identification numbers of the product are illegible or altered.

RESPONSIBILITY

The manufacturer or retailer decline any responsibility for any loss and direct or indirect damages, expenses that are directly or indirectly connected to the use of the appliance LCLC1330CE or because of consequences of events, such as diagnosis, medical procedures and treatments taken by the medical personnel that may cause damages to persons or objects, or cause expenses. Inside the packaging, we supply the form "Warning of Malfunctioning" that must be filled and sent to the manufacturer for any accident or near miss that may occur during the intervention.

The code LCLC1330CE is supplied with the connections ① and ② assembled and it is ready to be applied upon carriages. In case of use of the wall regulator, dismantle the connection ① (key 17) and assembly it again in place of the connection ② (key 14) and, in place of the connection ① assembly the desired wall connection ③.
Rev. 01 of 07/2019

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